


Run 8cm(h) x 10cm(w) SDS PAGE gels with unsurpassed resolution

FEATURES


- Triethanolamine (TEO) gel technology features a unique acrylamide polymer matrix for exceptional protein separations
- There are no combs or tape to remove shortening prep time.
- All gels have the ability to expand the resolving range of a given gel by running the small proteins off the unsealed bottom of the gel.
- Sharp, clean, consistently straight bands with unsurpassed resolution
- 40% faster separation with "Turbo" buffer
- Super-strong gel makes handling easy... no tearing, ripping or cracking. Ideal for Western blots.
- Simple buffer system. Reusable lower buffer saves time and money.
- TEO gels have a neutral pH for a 2 year shelf life - never throw away another gel
- Cost-effective replacement for Novex™, Invitrogen and Bio-Rad precast gels with a longer shelf life.

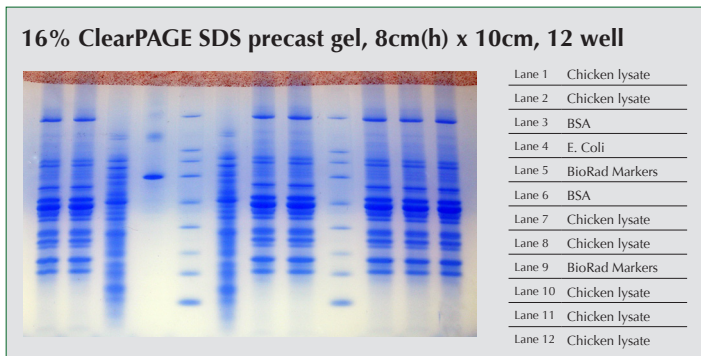


SPECIFICATIONS	
Gel types	Triethanolamine SDS gels
Shelf life	2 years
Overall cassette dimensions	8(h) x 10(w) x 5mm(d)
Gel dimensions	6.5cm(h) x 8.5cm(w) x 1mm(d)
Number of gels / box	10
Loading area color of TEO gels	red
Acrylamide %'s offered	8%
	10%
	12%
	4-8%
	4-12%
Number of wells offered	4-20%
	12 & 17

COMPATIBLE WITH:


C.B.S. Scientific Electrophoresis Systems

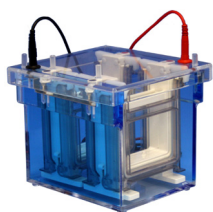

Quadra QNX-700
4-place Combo Cell

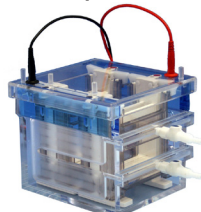


Quadra Cool QNC-700
4-place Combo Cell

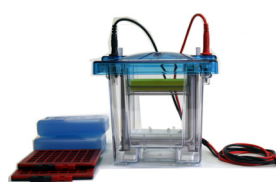

Dual Cool DCX-700
2-place Combo Cell
with freezer blocks



Dual Mini-Verticals
MGV-202 and MGV-
402, internal cooling

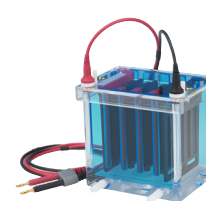

Mini-Blotter EBU-204
4-place, cooling base




Mini-Blotter EBX-700
4-place with freezer
blocks

Also compatible with:

- Lonza™ PAGER® Minigel Chamber
- Bio-Rad Mini-Protean® Systems (requires gasket # BRGK-100)

ClearPAGE™ is a trademark of C.B.S. Scientific. PAGER™ and Lonza™ are trademarks of Lonza Corporation. XCell™ Novex™ and Invitrogen™ are trademarks of Life Technology Corp. Mini-PROTEAN® is a registered trademark of Bio-Rad Laboratories Corporation.

8CM x 10CM CASSETTE FEATURES & DIMENSIONS

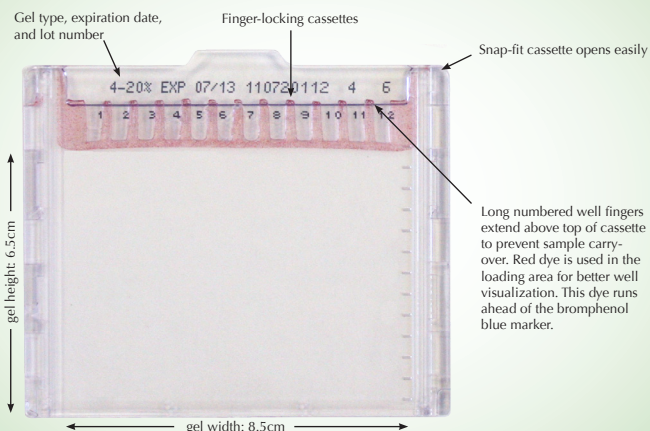

ORDERING INFORMATION


ClearPAGE™ SDS gels- TEO, 8cm(h) x 10cm(w)			
Cat. #	Quantity	Acrylamide %	# of wells
BK00812-10	Box of 10	8%	12
BK00817-10	Box of 10	8%	17
BK01012-10	Box of 10	10%	12
BK01017-10	Box of 10	10%	17
BK01212-10	Box of 10	12%	12
BK01217-10	Box of 10	12%	17
BK01612-10	Box of 10	16%	12
BK01617-10	Box of 10	16%	17
BK40812-10	Box of 10	4-8%	12
BK40817-10	Box of 10	4-8%	17
BK41212-10	Box of 10	4-12%	12
BK41217-10	Box of 10	4-12%	17
BK42012-10	Box of 10	4-20%	12
BK42017-10	Box of 10	4-20%	17

ClearPAGE™ buffers, stains and markers		
Cat. #	Description	Quantity
FB31002	LDS Sample Buffer	2 ml
FB31010	LDS Sample Buffer	10 ml
FB32001	DTT reducer, 10X	1 ml
FB50053	SDS-non-reducing running buffer	50 ml
FB50500	SDS-non-reducing running buffer	0.5L
FB60053	SDS reducing running buffer	50 ml
FB60500	SDS reducing running buffer	0.5 L
FB14053	Turbo buffer	50 ml
FB14050	Turbo buffer	0.5 L
FB82500	Transfer buffer	0.5 L
HG73010	Instant Blue stain	1 L
HG73050	Instant Blue stain, trial size	50 ml
HM05160	2-color marker	0.6 ml
HM05103	2-color marker, trial size	30 µl

ClearPAGE™ Blotting and Drying		
Cat. #	Description	Quantity
FB82500	TGS transfer buffer, 10x	0.5 L
HP19001	Blot sandwich NC 90 x 85mm, 0.2µl nitrocellulose/1mm blotting paper	2 pack
HP19020	Blot sandwich NC 90 x 85mm, 0.2µl nitrocellulose/1mm blotting paper	20 pack
HP29301	Blot sandwich PVDF 90 x 85mm, 0.2µl PVDF/1mm blotting paper	2 pack
HP29320	Blot sandwich PVDF 90 x 85mm, 0.2µl PVDF/1mm blotting paper	20 pack
HB04510	Gel drying solution	1 L
HB02520	Precut cellophane	20 sheets
HB90001	Drying frame with 8 clamps with cellophane	20 pieces of cellophane

Gasket Adapter for Mini-Protean® Systems		
Cat. #	Description	Quantity
BRGK-100	Gasket Adapter Kit for Mini-Protean®	1 pair

ClearPAGE™ 8cm(h) x 10cm(w) Precast Gel Migration Chart in kDa Turbo Buffer - Cat. # FB14500

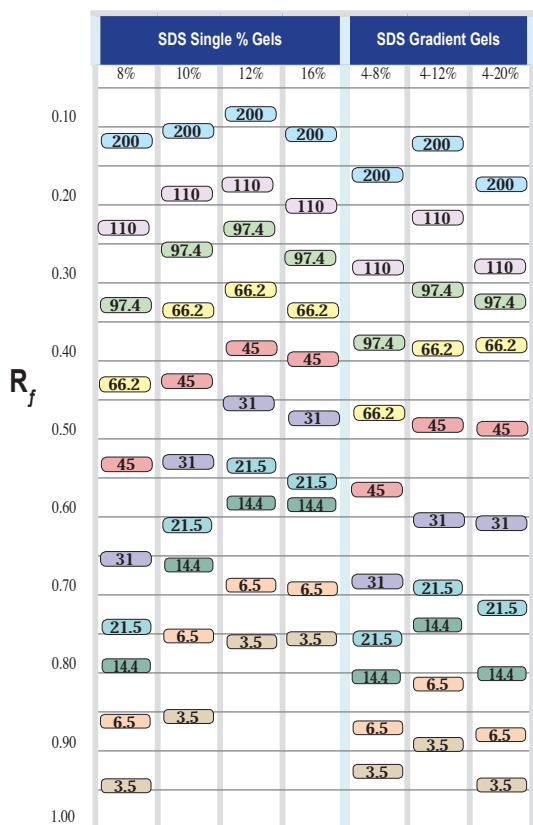

Gasket replacement is recommended for use of ClearPAGE™ 8 x 10cm Precast Gels in Mini-Protean® Systems

Easy to replace: simply remove the green gasket and replace with white gasket as shown at left.

FREE pair of gaskets with purchase of 1 box of gels if requested. Order gasket adapter kit BRGK-100